
Sleep.8

EN
PT
ES
LV
RO
EL
ET
FR
RU
DE
FI

POWER FOOT

Style, functionality and care for your feet!

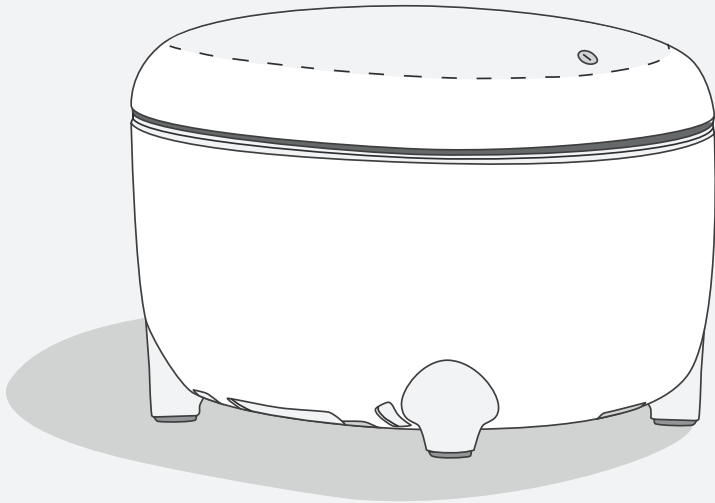Care and operations instruction


POWER FOOT

Style, functionality and care for your feet!

EN	2
PT	13
ES	24
LV	35
RO	46
EL	57
ET	68
FR	79
RU	90
DE	101
FI	112

POWER FOOT

Style, functionality and care for your feet!

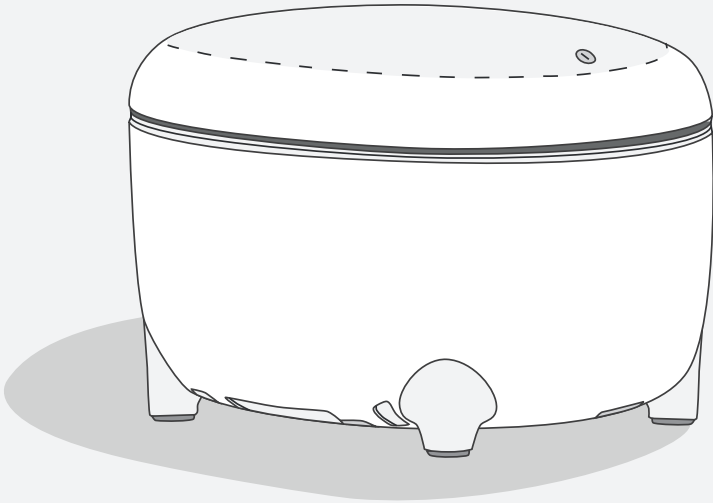Please read the Safety Precautions carefully before using the foot massager. The Safety Precautions are divided into ‘warnings’ and ‘attention’ in terms of the degree of danger. This book should be kept available for future reference after reading.

Safety Instruction	4
Product Parts	6
Product Features and Operation Illustration	7
Product Operation Guide	8
Product Maintenance	10
Troubleshooting	11
Product Specifications	12
Warranty	12

User precautions

- ▶ Make sure that this product is far away from any hot surfaces.
- ▶ Please consult your doctor before use for those who have the following diseases:
 - Those who have heart disease, malignant tumor, back pain, spinal scoliosis, hypertension, osteoporosis, etc.
 - Those who are pregnant or having a menstrual period.
 - Those who have a physical discomfort, fever, a swelling or inflammation, bruises, allergies, etc
 - Those who are warned by the doctor not to do any massage to avoid injury or physical discomforts such as blood clots (thrombosis), severe aneurysm, acute venous tumor, all kinds of skin inflammation, infectious skin disease (including subcutaneous tissue inflammation), etc.
 - Warnings: Do not use the massager together with the following medical instruments.
 - 1) Those who are using a pacemaker or other implantable medical equipment.
 - 2) Those who are using an artificial heart.
 - 3) Those who are wearing an electrocardiograph for similar medical equipment.

Otherwise it will cause malfunction for those medical equipment, putting the user's life in danger.

Precautions for power supply

- ▶ Do not turn on the switch when your hands are wet.
- ▶ Do not immerse in or spray water or liquid on the control power switch.

Cord/Plug Precautions

- ▶ Do not connect or remove the plug with wet hands. Otherwise you will suffer an electric shock or a product malfunction.
- ▶ Do not pull the wires when you pull out the power plug. Please hold the plug itself.
- ▶ As the default operating time is 20 minutes, all the functions will stop working after 15 minutes. Please unplug this product from the electrical outlet immediately after use.
- ▶ Keep the cord away from hot surfaces.
- ▶ Do not place the power plug under the product or any heavy items.

SAFETY INSTRUCTIONS


- ▶ Do not wind the power plug around the product lest it should be damaged and cause an electric shock or fire.
- ▶ Do not operate this product if it has a damaged cord or plug.
- ▶ Do not use this product if the power socket is loose.
- ▶ Do not carry this product by its supply cord or use the cord as a handle.
- ▶ When the power plug is damaged, please call the service centre or an authorized technical person for the necessary check and maintenance.

Attention



- ▶ Do not fall asleep while massager is on.
- ▶ Do not use this product while you are drunk or physically uncomfortable.
- ▶ Don't stand on the product or put a heavy weight on it.
- ▶ Do not use this product for more than 30 minutes.
- ▶ Do not use this product if the cloth cover is worn out.
- ▶ Do not use this product within an hour before or after a meal.
- ▶ When you are using the product at a low room temperature, do not suddenly use the heat function. You can adjust the room temperature to normal room temperature, before starting the heat function.
- ▶ When moving this product from a low temperature storage place to a warmer setting, it may cause malfunction due to water vapor gathered in the components. Set the product in a normal room temperature for one hour before using it.
- ▶ Do not use this product continuously if you find any components damaged or the product in an abnormal condition.
- ▶ When current interruption occurs during your massage, try to avoid further damage due to sudden power restoration.
- ▶ Stop using the product if there is thunder and lightning.

PRODUCT PARTS

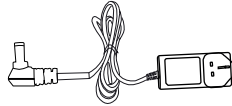The power adapter

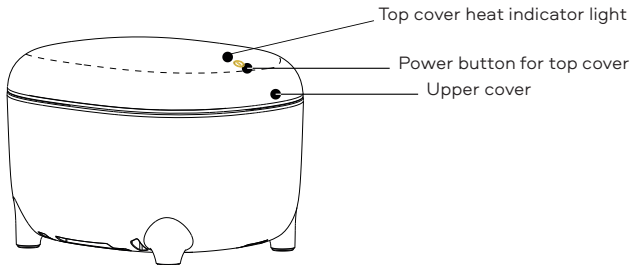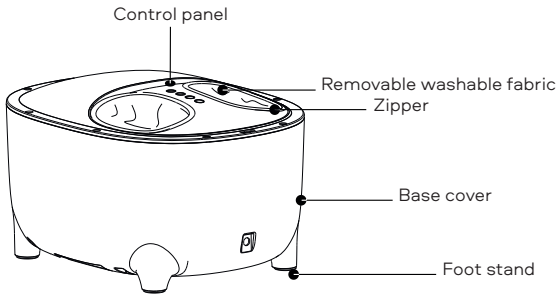
Preparation before use

- ▶ Please put the product on the at floor, and then remove top cover.

Usage method

Step 1 : wire connection.

- ▶ 1. Connect power adaptor to the power line plug-hole located on the product as shown on the right.
- ▶ 2. Plug the power plug into the AC 100-240~V power outlet.

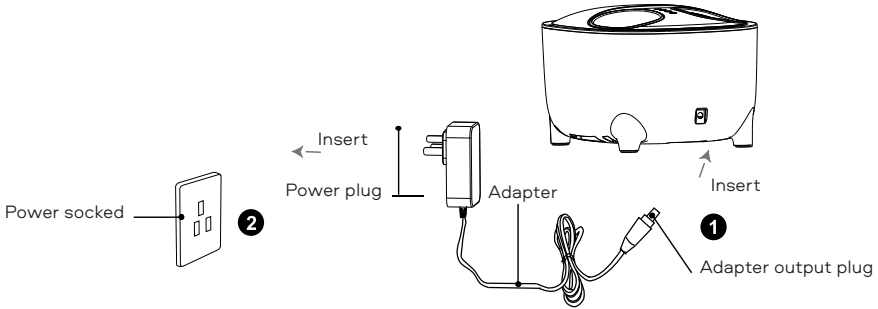⚠ CAUTION

- ▶ Before plug, please check make sure the power cord is in good condition.
- ▶ Do not use the wet hand to pull out the power cord inorder to avoid any electric shock.

Step 2 :

- ▶ Please choose a chair with desired height and then put your feet into the product (Picture 1)

Step 3 :

- ▶ Please turn on the power button «»
- ▶ The upper cover can be used as a stool (Picture 2). If want to be heated, press the heating power button on the cover, which will begin to heat (Picture 3).

Picture 1

Picture 2

Picture 3

After use

- ▶ Turn of the power button on the control panel.
- ▶ Remove the plug from the power supply socket when not in use.

⚠ CAUTION

- ▶ Please do not touch power connect points after cover heated in case of a hot.

ON/OFF

- ▶ Press to turn ON.
- ▶ Press again to turn OFF.
- ▶ Default setting timer is 15 minutes.

Auto mode button

Press to enable AUTO program 1 (green light)

- ▶ foot air bag+heat
Press again to switch to program 2 (blue light)
- ▶ foot roller+heat+foot tapping
Press again to switch to program 3 (red light)
- ▶ foot air bag+foot roller+heat+foot tapping
Press again to switch to program 4 (purple light)
- ▶ foot air bag+foot roller+heat+foot tapping+ankle air bag

Heat button

Press to turn ON.
Press again to turn OFF.

Intensity button

Press to switch message INTENSITY.
Weak (green) > Middle (blue) > Strong (red)

Power button on top cover

Press to turn ON heating function on top cover.
Press again to turn OFF.

CAUTION

- ▶ Please do not use the product by standing posture or put any heavy object on the product to avoid any damage of the product.
Please do not stand or seat on the product.
- ▶ When moving the product, please do not pull the power cord or drag the product directly to avoid any damage.
- ▶ Do not try the product on a no-load basis to avoid wear and tear.

- ▶ Each massage program is 15 minutes. Product will stop once massage is completed.


Instructions before use

Please store this product in a dry place, do not expose the product to sunlight put the carpet under the product to minimize the working noise.

How to use it

This product uses only the AC 220V-240V, 50/60Hz power supply. Make sure your indoor voltage fits this specification parameter.

- ▶ Connect the AC power cord to the socket.

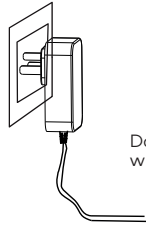Do not pull the plug
with wet hands.


Storage

- ▶ Please store it in the packing carton, and put in a safe, dry, well-ventilated, and shaded place.
- ▶ Do not twist the power plug.
- ▶ Do not hand the power wire and pull the product.
- ▶ Avoid sharp items from destroying the surface of the product.
- ▶ Please don't store the machine under sunlight or in places of high temperature.

Clean

- ▶ Before cleaning, make sure that the power cord is unplugged from the electrical outlet and that the machine is completely cool.
- ▶ Use a damp cloth to clean the product. Ensure the product is completely dry before reuse.
- ▶ The detachable fabric cloth can be washed in a normal washing way.

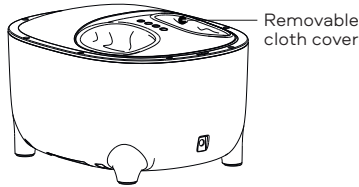Attention

- ▶ Make sure not to use benzene or thinner and other corrosive detergents to clean the product.

Maintenance

- ▶ Keep the product clean and store in a dry environment.
- ▶ When the product is not in use for a long time, please put a dustproof cover on the product.

If you meet the following problems, please referring to the following solutions. If the product is still malfunctioning or doesn't work, please send the product to our service centre or an authorized technician for the necessary checks and repair.

Conditions	Possible reason	Solutions
It doesn't work	<ol style="list-style-type: none">1. Plug isn't connected correctly to the socket.2. The power switch at the back is off.3. Switch key on control panel isn't on.	<ol style="list-style-type: none">1. Push the plug firmly into the socket.2. Switch on the power switch at the back of the product.3. Switch on the on/off key on the panel.
Stop working suddenly	<ol style="list-style-type: none">1. The power switch is off.2. Massage time is over 15 minutes.	<ol style="list-style-type: none">1. Check that the power plug is connected well.2. Restart the power supply switch.

PRODUCT SPECIFICATION


Power supply	AC 220V-240V, 50/60HZ
Power consumption	50W
Timer	15 minutes
Dimensions (LxWxH)	454.5x358.5x300 mm
Weight	7.3 Kg
Dimensions of packing box (LxWxH)	485x338x388 mm
Weight (including the packing box)	8 Kg

WARRANTY

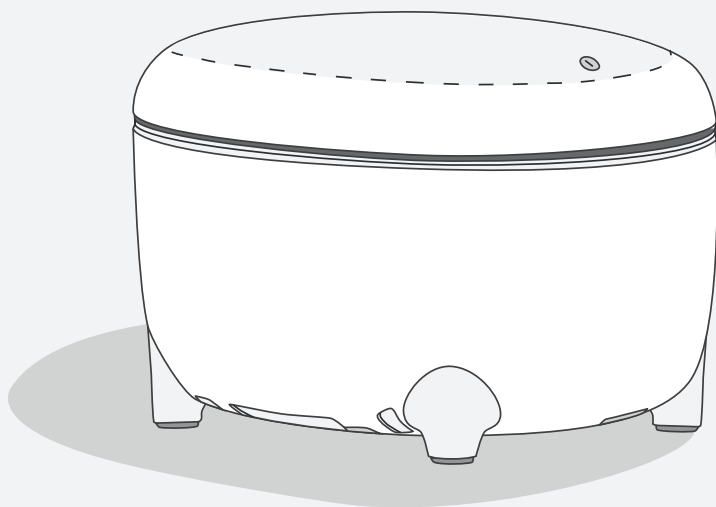1. The product is guaranteed for 12 months from the date of purchase.
2. The warranty is valid if the product is defective due to the fault of the manufacturer.
3. The warranty does not apply to mechanical damage and if there are signs of repair by unauthorized persons.

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