
EN
PT
ES
LV
RO
EL
ET
FR
HE
DE
RU

STRESS FREE

*Indulgent facial massage from
the comfort of your home*

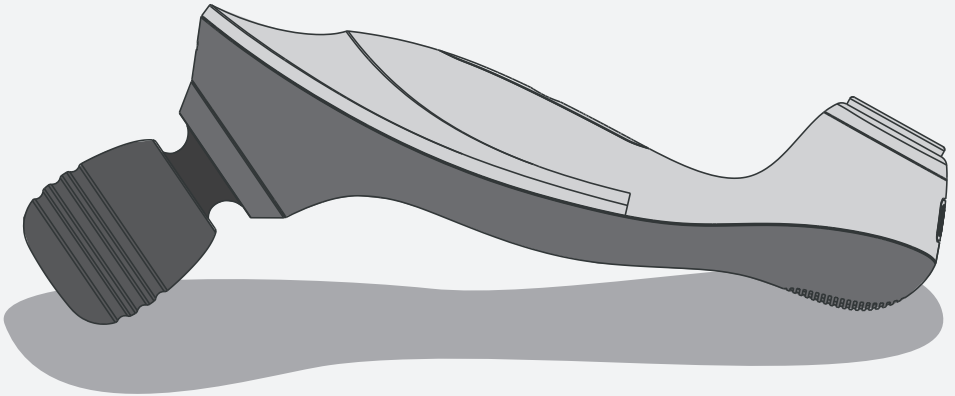STRESS FREE


*Indulgent facial massage from
the comfort of your home*

ENGLISH	2
PORTUGUESE	6
ESPAÑOL	10
LATVIAN	14
ROMANIAN	18
GREEK	22
ESTONIAN	26
FRENCH	30
HEBREW	34
GERMAN	38
RUSSIAN	42

SAFETY PRECAUTIONS

1. The following conditions should consult doctor before use:
 - people suffer from heart disease
 - people suffer from infectious disease
 - people suffer from malignant tumor
 - people has open cut on the skin
 - women in pregnancy, menstruation and perpetual fever
 - people with excessive sunburn
 - people with sensitive skin or angiectasis or have been admitted to hospital.
2. Do not disassemble the device by yourself. Seek help from the authorized service center when necessary.
3. Do not use or place the device in a hot or humid environment to prevent functional defect or damage.
4. Do not charge the device under wet area to prevent the short-circuits.
5. Do not allow children or people with physical limitations to use the device without supervision.
6. This device is not water-proof. Do not leave it under water or inside corrosive environment.
7. Stop using the device immediately when any abnormalities from the device occurred.
8. Do not drop or hit the device to avoid any damage.
9. Essential oils may be applied during using the devlce to prevent skin from getting too dry.

PRODUCT DESCRIPTION

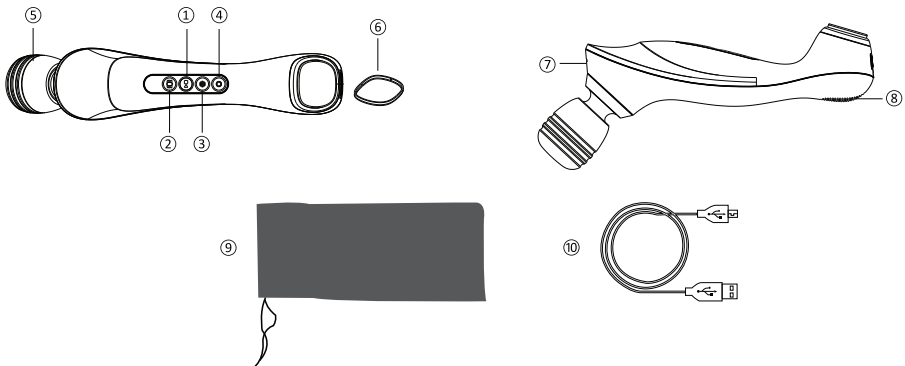- | | |
|-----------------------|--|
| 1. Main power switch | 6. Protective Cover (hot/cold) alternating metal sheet |
| 2. Intensity adjuster | 7. USB charging port |
| 3. Cool Mode | 8. Convection window |
| 4. Warm Mode | 9. Storage bag |
| 5. Massage Head | 10. USB Cable |

HOW TO OPERATE

Before Use. Please fully charge the device before using it for the first time.

HOW TO CHARGE

1. Connect one end of the USB cord provided to the charging port and insert the other end to the external power source.
2. While charging the LED light beneath the power button will turn red; when fully charged the LED light will turn violet. Estimated charging time is between 4-6 hours.
3. During charging process, the device will not be able to function as it has a self-protection device installed.
4. This device uses high-performance & environment-friendly lithium ion battery; the life of the battery varies according to the function of the product. Please refer to the manual for details.

When the battery power is low, the LED light beneath the power button  will blink thereafter the device will shut down automatically.

VIBRATING MESSAGE

1. Press the power switch  and the vibration will begin
2. 3 massage intensities are available for selection; press and hold the vibration button  the intensity level will switch accordingly.
3. The default time for vibration massage is 10 minutes; press the power button again to continue the massage.

COLD THERAPY

Remove the protective cover on the contact metal of the device.

1. Press the power button  to activate the device.
2. Press the cold therapy button  to activate the cold therapy mode
3. There are three levels of cold temperature for selection; press and hold the cold therapy button  the temperature levels will switch accordingly.
4. The temperature levels are indicated by the following LED light beneath the cold therapy button:
 - a. Blue represents level 1 (approximately 15-20°C)
 - b. Violet represents level 2 (approximately 10-15°C)
 - c. Red represents level 3 (approximately 5-10°C)
5. The default time for the cold therapy is 3 minutes, to continue press the cold therapy button  again.

HEAT THERAPY

Remove the protective cover on the contact metal of the device.

1. Press the power button  to activate the device.
2. There are three levels of heat temperature for selection; press and hold the heat therapy button  to switch accordingly.
3. The temperature levels are indicated by the following LED light beneath the heat therapy button:
 - a. Blue represents level 1 (approximately 30-35°C)
 - b. Violet represents level 2 (approximately 35-40°C)
 - c. Red represents level 3 (approximately 40-45°C)
4. The default time for heat therapy is 3 minutes; to continue press the heat therapy button  again.

CARE & MAINTENANCE

- Always unplug from power cable and switch off after use.
- Put back the protective cover onto the metal contact.
- Clean the device with a soft & dry cloth when necessary.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, glass furniture polish or thinner on the device.
- Store the device in a dry, cool place.
- Do not put heavy objects onto the device.

TROUBLESHOOTING

This page summarizes the common problems you may encounter while using the device. If you are unable to resolve the problems with the suggestions below, please contact the authorized service center.

Problem	Possible Cause	Suggestion
Product not functioning	Product not sufficiently charged	Charge the product
	Default-time is up	Restart
Product is not charging	Power supply not connected correctly	Unplug and re-connect
	Not using the USB cable provided	Restart
Not fully charged even after long charging time (the indicator stays on red)	Power source has insufficient voltage	Change the power source or charge from computer or power bank. Different charging source required different charging time.
Slight sound heard during operation	Sound is the result of mechanical movement	Normal phenomenon. Product not faulty. No action required

PRODUCT SPECIFICATIONS

Product Name:	Cold & Hot Beauty Massager
Input Voltage:	USB 5V, 1000mAh
Battery parameters:	3.7V 2600mAh
Temperature range:	5-45°C (±5°C)
Power MAX:	10W
Charging time:	4 - 6 hours
Motor Speed:	4000-8000 RPM
Battery life:	<ul style="list-style-type: none">◦ 1,5 - 2,5 hours with vibration on◦ 1,5 - 2,5 hours with cooling◦ 4 - 5 hours with heating
Automatic Timer:	<ul style="list-style-type: none">◦ 10 min with vibration◦ 3 min with cooling◦ 3 min with heating
Product Dimensions:	280x65x60 mm
Net weight:	370 g
Safety Compliance:	CE EN60335-1, EN60335-2-32

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