
EN  
PT  
ES  
LV  
RO  
EL  
ET  
FR  
HE  
DE  
RU

## POWER EYES

*Like a spa day for your eyes*
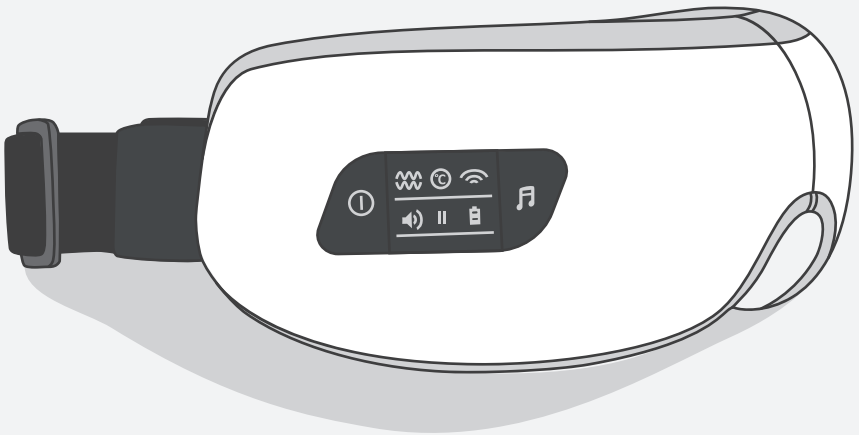
# POWER EYES

*Like a spa day for your eyes*

ENGLISH	2
PORTUGUESE	7
ESPAÑOL	12
LATVIAN	17
ROMANIAN	22
GREEK	27
ESTONIAN	32
FRENCH	37
HEBREW	42
GERMAN	47
RUSSIAN	52

## Power Eyes Massager


Sleep8 Power Eyes Massager is a solution specially designed to relieve tension of your eyes and also for rest and recreation after a busy day. The wireless gadget can be used at home or on the move. You can take it with you in a bag due to its compact design.

### THE ADVANTAGES OF THIS MASSAGER

- It stimulates blood circulation.
- It helps to relieve headache.
- It relieves muscle tension while improving flexibility.
- It shortens the restoration period.


Please, read the current manual carefully before using the Massager.

## Description of Power Eyes

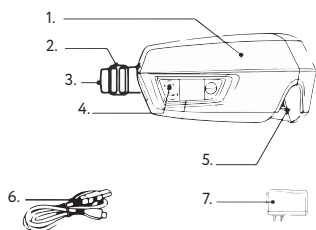1. Case
2. Size Plate
3. Elastic Band
4. Control Panel
5. Air Pad
6. USB Cable for Charging and Music Recording
7. Adapter

## Control Panel

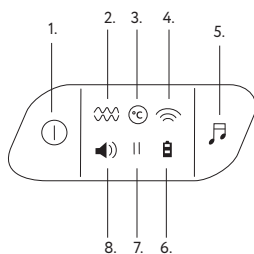1. Device On Button
2. Air Compressor Indicator
3. Warming Mode Indicator
4. Vibromassage Multisync Indicator
5. Music Control Button
6. Battery Indicator
7. Playing Music Indicator
8. Music Playback Indicator

## Safety Precautions

1. Before using, please consult a doctor if:

- You suffer from eye diseases such as glaucoma, leukoma and retinal detachment.
- You suffer from various weaknesses and fatigues, and also physical, mental or psychogenic disorders.
- You suffer from insult, hypertonia/hypotonia or heart diseases.
- You suffer from infectious diseases.
- You suffer from a cancer.
- You have had an eye surgery.
- You have a diagnosis of schizophrenia.
- You have open wounds, cuts, redness, infections, burns or bruises on the skin.
- You are pregnant, in the period of menstruation and incessant fever.
- You have a sensitive skin, angiectasia or were hospitalized.
- you have low temperature sensitivity, maximum temperature of the warming mode is 42°C.

2. Do not use the Massager when driving a car.

3. Do not disassemble the device yourself. If necessary, contact the service center.

4. Do not use, charge or store the device in hot or humid environment or in direct sunlight in order to avoid failure and malfunction.

5. Do not allow children and disabled people to use the device unattended.

6. Do not leave the massager under water or in rooms that provoke corrosion.

7. In case of malfunction stop using the device immediately.

8. Do not drop or hit the device.

9. Do not put heavy objects on the massager.

10. The massager is not a replacement for treatment or doctor's recommendations. If you are in a doubt please consult your doctor.

11. The massager shall not be used for diagnosis, treatment or for commercial purposes.

12. Do not use the massager during charging.

# Massager Operating Instructions

## BEFORE USE

- Check the availability of the device and components.
- Fully charge the massager.
- Be sure to take out contact lenses and eyeglasses.

Turn on the massager only after you put it on your eyes in order not to damage the air pad due to over-pressure.

Do not start the device if the power and power supply voltage are out of range or stay unstable.

## HOW TO CHARGE THE DEVICE

1. Connect the USB port of the massager to an external power source using the USB cable.
2. Connect the charger to the socket (100-240 V ~, 50/60 Hz). When the device is connected to the charger the LED will turn blue.
3. The charge will take 3 hours.

## OPERATING PROCEDURES

1. Put the Massager on your eyes.
2. Be sure that the Massager fits tightly to your face but does not cause any discomfort. You can adjust the way it fits by using a buckle on the elastic strap.
3. Start the device by pressing and holding the button  : After the switch on the LED panel goes purple. The massage with warming effect and playing music will start automatically after switch on.
4. To complete using the Massager press and hold the  button. The LED blinks off.
5. The Massager automatically stops working in 15 minutes. To resume work press and hold the  button.

## POSSIBLE MESSAGE MODE RANGE

1. Delicate air-compression massage.
2. Vibromassage of moderate intensity.
3. All-in massage with interleaving air-compression and vibromassage supplemented by the light warming effect.
4. The default mode is air-compression massage with the light warming effect.

# Massager Operating Instructions

## HOW TO CHOOSE THE DESIRED MODE

1. The default mode starts by pressing and holding the button switch on when the Massager is switched off.
2. The vibromassage mode starts by one quick pressing when the default mode is running.
3. The all-in mode starts by quickly pressing the button twice when the default mode is running.
4. Air-compression massage starts by quickly pressing the button three times when the default mode is running.

## PLAYING MUSIC FUNCTION

To start the Playing Music function do the following:

1. Press and hold the Play button  .
2. In order to switch a song, press the button to play a song twice.

## Troubleshooting

When using the Massager, you can experience the following. If these recommendations have no desired effect or are not applicable please contact the service center.

Type	Possible Reason	Troubleshooting
The device is not working	Check if the device is switched on	Press the ON/OFF button
	Low battery charge	Please charge the Massager
The device fails to charge	The power supply is not connected properly	Disconnect and then connect the Massager to the power supply
The LED is flashing	Low battery charge	Please charge the Massager

## **Clean and Care**

### AFTER USAGE

- Turn off the Massager and disconnect the USB cable.
- Store in a dry cool place. DO NOT store the device up in the air.

### CLEAN

- Wipe the Massager with a soft dry cloth. Do not use water to clean the device.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, furniture polishes or their solutions.

## **Technical Specifications of the Product**

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