
EN
PT
ES
LV
RO
EL
ET
FR
HE
DE
RU

POWER NECK

Deep, soothing relaxation for your neck

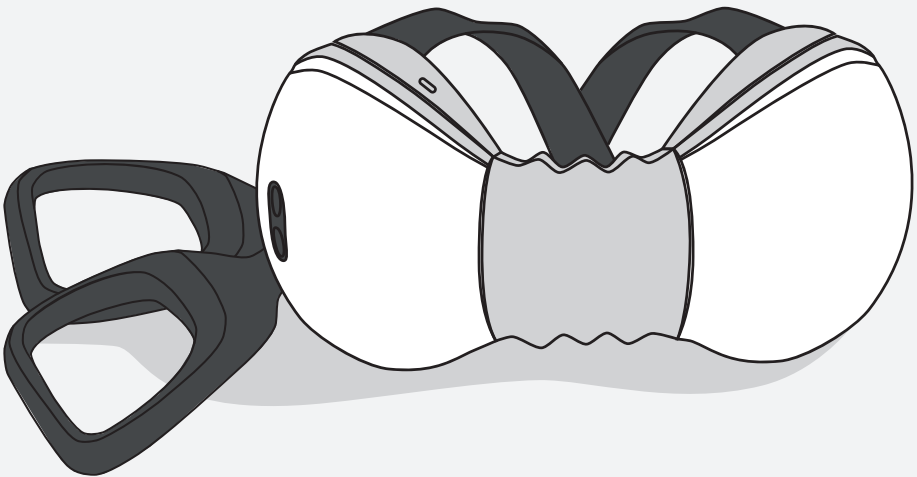POWER NECK

Deep, soothing relaxation for your neck

ENGLISH	2
PORTUGUESE	7
ESPAÑOL	12
LATVIAN	17
ROMANIAN	22
GREEK	27
ESTONIAN	32
FRENCH	37
HEBREW	42
GERMAN	47
RUSSIAN	52

POWER NECK MASSAGER


Power Neck Massager is a gadget specially designed to relieve tension of neck and collar zone muscles. Sleep 8 Power Neck reproduces a traditional Japanese massage Shiatsu technique. In addition to the mechanical effect of the massager it has a light warming mode to improve impact on muscles. After using this device, you will feel that your mood improves, your productivity increases and the tension of your muscles relieves.

THIS DEVICE IS DESIGNED SPECIALLY FOR THOSE WHO:

1. Spends more than a couple of hours per day in front of a computer and stays in one position for a long period of time.
2. Suffers from nervous tension, insomnia or sleepiness.
3. Suffers from neck and shoulder tension and discomfort.
4. Suffers from a lack of oxygen supply to the brain, headache and faintness.

Before using, please read this User Manual.

DESCRIPTION OF POWER NECK

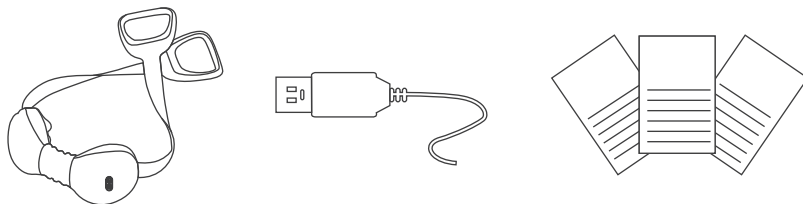-
1. Device On Button
 2. Warming Mode On Button
 3. A Massage Head
 4. A USB Charging Port
 5. A USB Cable

SAFETY PRECAUTIONS

1. Before using, please consult a doctor if:

- you suffer from heart diseases or high blood pressure
- you suffer from infectious diseases
- you suffer from a cancer
- you have a diagnosis of schizophrenia
- you have an open wound, cut, redness, infections, burn or bruise on the skin
- you are pregnant, in the period of menstruation and incessant fever
- you have a sensitive skin, angiectasia or were hospitalized
- you had a neck surgery
- your joints and muscles can be described as slow to warm up

2. Do not use the Massager when driving a car.

3. Do not disassemble the device yourself. If necessary, contact the service center.

4. Do not use, charge or store the device in hot or humid environment or in direct sunlight in order to avoid failure and malfunction.

5. Do not allow children and disabled people to use the device unattended.

6. Do not leave the massager under water or in rooms that provoke corrosion.

7. In case of malfunction stop using the device immediately.

8. Do not drop or hit the device.

9. Do not put heavy objects on the massager.

10. The massager is not a replacement for treatment or doctor's recommendations.

11. The massager shall not be used for diagnosis, treatment or for commercial purposes.

12. Do not use the massager during charging.

MASSAGER OPERATING INSTRUCTIONS

- Check the availability of the device and components.
- Fully charge the massager before use.
- If you have long hair gather it on top of your head before use.
- Turn on the massager only after you put it on your neck.
- Do not start the device if the power and power supply voltage are out of range.

HOW TO CHARGE THE DEVICE

1. Connect the USB port of the massager to an external power source using the USB cable.
2. When charging, the LED will turn red. When the device is fully charged this indicator will turn purple. The device uses a highly effective lithium battery.
3. The battery is a highly effective lithium battery.

OPERATING PROCEDURES

1. Hold the handles with both hands, put the Massager on the cervical spine, adjust and find convenient position suitable for your neck.
2. Pull the handles with both hands and adjust position of the massage heads.
3. Start the device by pressing and holding the button  : After the device is switched on the LED goes purple.
4. To change direction, press the  button.
5. To complete usage press and hold the  button. The LED blinks off.
6. Heating is enabled automatically when the massager is turned on. The  button turns the heating on and off.
7. The massager automatically stops working in 20 minutes. To resume work press and hold the  button.
8. Flashing red LED informs about low battery charge. Please charge the device.

TROUBLESHOOTING

When using the massager, you can experience the following. If these recommendations have no desired effect or are not applicable please contact the service center.

FAILURE	POSSIBLE REASON	TROUBLESHOOTING
The massager does not switch on	Low battery charge	Charge the device fully
	Incorrect start of the Massager	Press and hold ON/OFF to turn on the device
The device fails to charge in case of long charging (the LED indicator of the battery is red)	The voltage of the USB charger may not correspond to the main device	Replace the charger (mobile charging unit of any brand) or connect to the computer, power supply unit
The device stopped working during massage	Low battery charge	Please charge the massager
	LED indicator is red, flashes 8 times and stops the device	Restart the massager
Slight massage effect	Check the battery level	Charge the device

CLEAN AND CARE

AFTER USAGE:

- Turn off the massager and disconnect the USB cable.
- Store in a dry cool place. DO NOT store the device up in the air.

CLEAN:

- Wipe the massager with a soft dry cloth. Do not use water to clean the device.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, furniture polishes or their solutions.

TECHNICAL SPECIFICATIONS OF THE PRODUCT

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